

INTERVENTIONAL
VASCULAR
DIAGNOSTICS &
THERAPY

VENA CAVA FILTER

VenaTech® Convertible, VenaTech® LP, VenaTech® Retrievable

B. BRAUN VENA CAVA FILTERS

For more than 20 years, inferior vena cava filters (IVC) have provided a protection of the respiratory function by decreasing the risk of pulmonary embolism.


B. Braun offers a complete range of permanent , retrievable and convertible IVCs.

ALL FILTERS HAVE A PROVEN:


- Clinical design for enhanced clot trapping*
- Long stabilizing filter legs, e. g. VenaTech Convertible*
- Combined delivery system for jugular and femoral approach
- Brachial approach available for permanent and retrievable filters

PRODUCT FEATURES:

- Made of non-magnetic cobalt-chromium alloy
- X-ray visible
- MR conditional

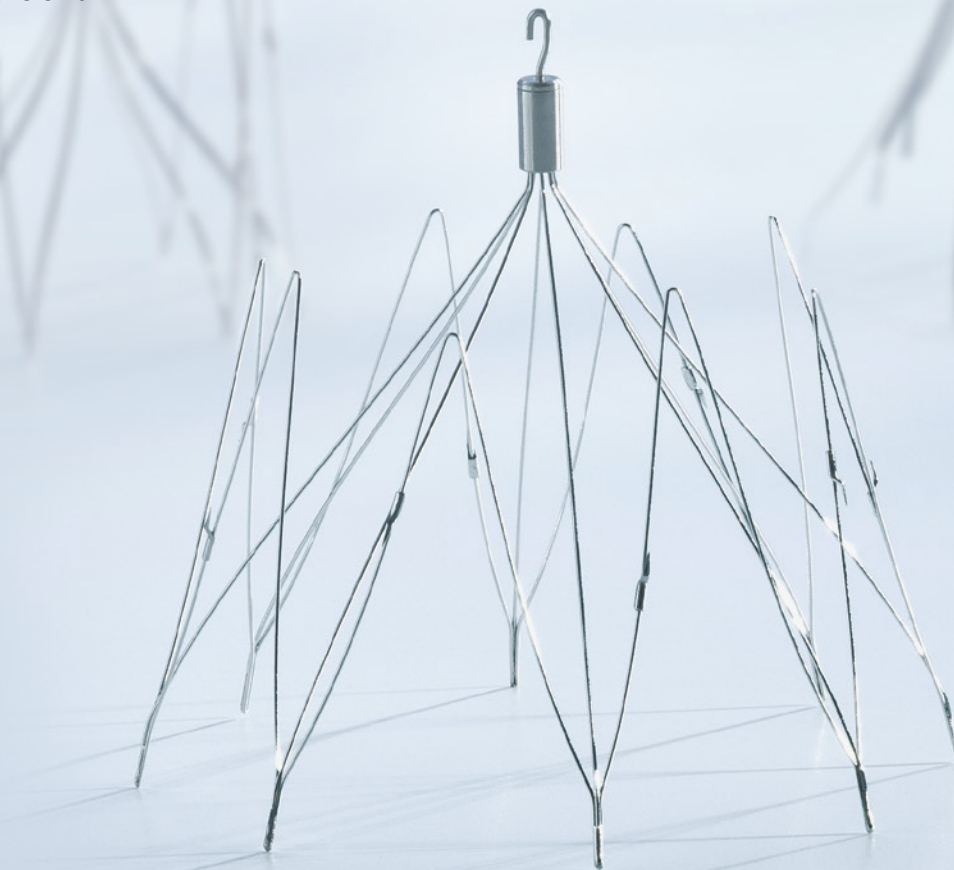


CONVERTIBLE VENA CAVA FILTER

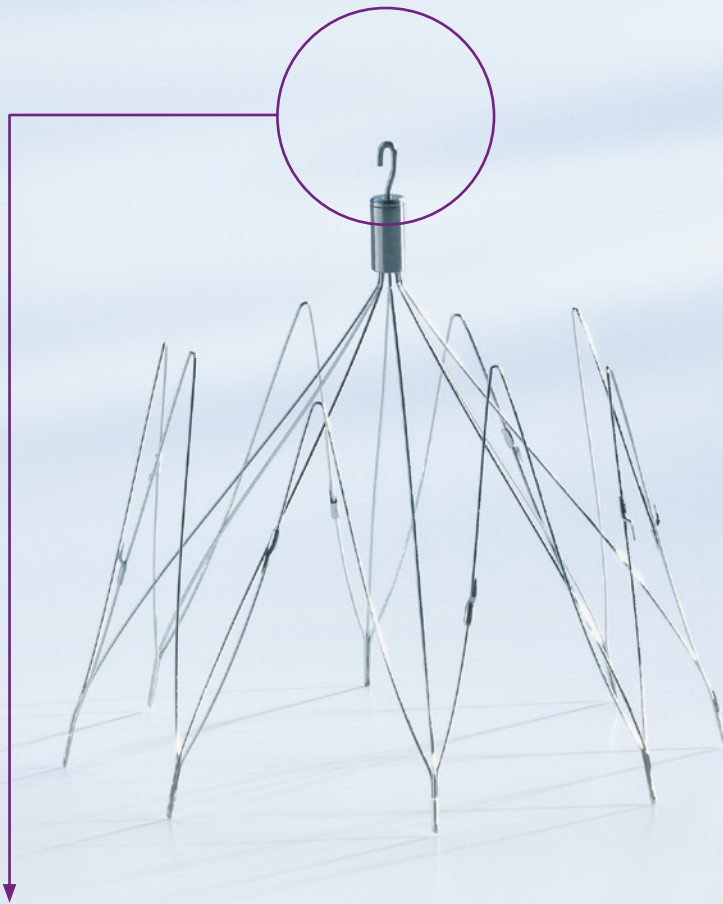


Filter Deactivation

VenaTech® Convertible - The Convertible Vena Cava Filter


- Effectiveness of caval filtration is combined with the possibility of filter deactivation.
- For femoral or jugular approach.
- Wire design and high flexibility to accommodate vena cava up to 32 mm in diameter.
- Patented concept of deactivation:
In cases where the physician determines that it is clinically indicated to discontinue filtration, the intraluminal filter elements of the VenaTech® Convertible vena cava filter can be converted to an open configuration. The time period to conversion for clinical study subjects ranged from 15 days to 391 days.

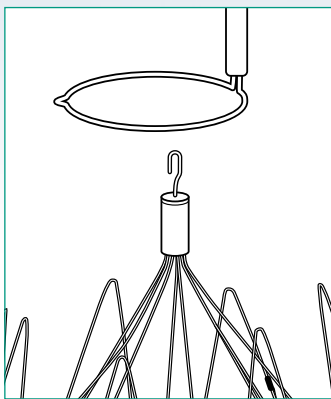


1 Advance gooseneck snare catheter to the top of the filter hook

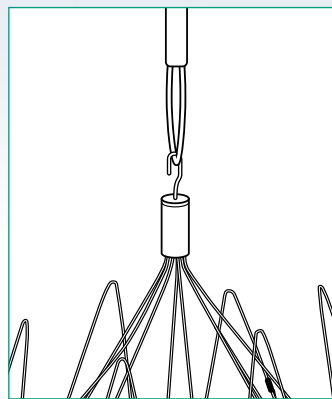

2 Capture the hook on the cone of the filter with the loop of the gooseneck snare

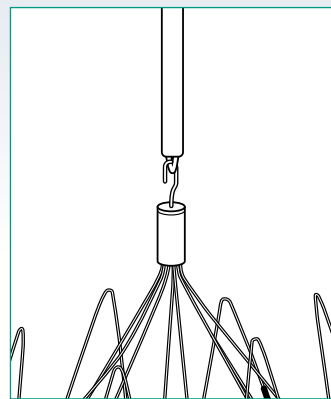

3 Advance the gooseneck catheter downward until it covers the hook

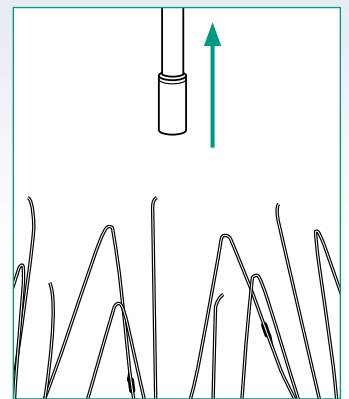

4 Retrieve filter head

- Pin the catheter in place and pull the snare proximally until the filter head unlocks.

PERMANENT AND TEMPORARY VENA CAVA FILTERS



VenaTech® LP

VenaTech® LP Permanent Vena Cava Filter

Features and Advantages

Long-term protection against recurrent pulmonary embolism

Improved implantation procedure

- For femoral, jugular, subclavian or brachial approach, you always have the right system available.
- Each filter cartridge is uniquely marked and designed to prevent misloading of the filter into the introducer sheath.
- VenaTech® LP has self centering, stabilizing legs which avoid filter tilting during the implantation.

Enhanced safety & stability

VenaTech® LP has patented stabilizing legs with securing hooks which were designed and calibrated to reduce the risk of migration while maintaining the integrity of the vena cava wall.



VenaTech® Retrievable

VenaTech® Retrievable Temporary Vena Cava Filter System

For decades retrievable vena cava filters have been well established in the market. B. Braun analyzed the existing retrievable filter concepts and developed a new filter design. Based on an established eight leg filter architecture, VenaTech® Retrievable was designed to have four horizontal levels of vessel contact to minimise tilting (4-level-safety concept).

Implantation period

The VenaTech® Retrievable, Vena Cava Filter System can be safely retrieved up to 12 weeks* after implantation. After this period of time, retrieval may be impossible and the permanent indication must be considered.

Four good reasons to use a VenaTech® Retrievable:

- Improved tilt resistance
- Safety cartridge to reduce the risk of sharps injury during filter retrieval
- Reduced penetration risk
- Ease of retrieval

* Janjua M, Younas F, Moinuddin I, et al. Outcomes with retrievable inferior vena cava filters. J Invasive Cardiol. 2010;22(5):235-239. PMID: 20440042

Brown JD, Raissi D, Han Q, Adams VR, Talbert JC. Vena Cava Filter Retrieval Rates and Factors Associated With Retrieval in a Large US Cohort. J Am Heart Assoc. 2017;6(9):e006708. Published 2017 Sep 4. doi:10.1161/JAHA.117.006708 PMID: 28871041

ACR-SIR-SPR PRACTICE PARAMETER FOR THE PERFORMANCE OF INFERIOR VENA CAVA (IVC) FILTER PLACEMENT FOR THE PREVENTION OF PULMONARY EMBOLISM - Revised 2021 (Resolution 8)

Clinical Study Documentation of the VenaTech® Retrieval System, Vena Cava Filter, Clinical Report, B. Braun Medical SAS

ORDERING INFORMATION

VenaTech® Convertible

Filter	Components	Reference
VenaTech® Convertible Convertible Vena Cava Filter System Femoral / Jugular	<ul style="list-style-type: none"> ▪ Introducer sheath (13 F O. D.*) ▪ Pre-loaded filter ▪ Pusher (10 F) ▪ Dilator (10 F) ▪ PTFE coated "J" guidewire, 1500 mm, 0.035" (0.89 mm) 	4435140

VenaTech® LP

Filter	Components	Reference
VenaTech® LP Permanent Vena Cava Filter System Jugular / Femoral	<ul style="list-style-type: none"> ▪ Vena Cava diameter < 35 mm ▪ Introducer sheath (7 F I. D.*/9 F O. D.*), hemostasis valve, 560 mm usable length ▪ Pre-loaded filter ▪ PTFE coated "J" guidewire, 1500 mm, 0.035" (0.89 mm) ▪ Long and short pushers (7 F) ▪ Dilator (7 F) 	4435125
VenaTech® LP Brachial Introducer System Introducer sheath (7 F I. D./9 F O. D.) (Antecubital)	<ul style="list-style-type: none"> ▪ Introducer sheath (7 F I. D.*/9 F O. D.*), hemostasis valve, 960 mm usable length ▪ Long and short pushers (7 F) ▪ Dilator (7 F) 	4439985

VenaTech® Retrievable

Filter	Components	Reference
VenaTech® Retrievable Temporary Vena Cava Filter System Jugular/Femoral	<ul style="list-style-type: none"> ▪ Vena Cava diameter >14mm and < 28 mm ▪ Introducer sheath (7 F I. D.*/9 F O. D.*) hemostasis valve, 560 mm usable length ▪ Pre-loaded filter ▪ PTFE coated "J" guidewire, 1500 mm, 0.035" (0.89 mm) ▪ Long and short pushers (7 F) ▪ Dilator (7 F) 	4435150
VenaTech® Retrievable Brachial Introducer System (Antecubial)	<ul style="list-style-type: none"> ▪ Introducer sheath (7 F I. D.* 9 F O. D.*) hemostasis valve, 960 mm (usable length) ▪ PTFE coated "J" guidewire, 1800 mm, 0.035" (0.89 mm) ▪ Long and short pushers (7 F) ▪ Dilator (7 F) 	4435160
VenaTech® Retrievable Vena Cava Filter Retrieval System Jugular	<ul style="list-style-type: none"> ▪ Introducer sheath (13 F O. D.*) ▪ Retrieval cartridge ▪ PTFE coated "J" guidewire, 1500 mm, 0.035" (0.89 mm) ▪ Dilator (7 F) 	4435170

* O.D. = outer diameter, I.D. = inner diameter

B. Braun Melsungen AG | Vascular Systems | Sieversufer 8 | 12359 Berlin | Germany
Phone +49 30 568207-300 | Fax +49 30 568207-210 | www.bbraun.com

Information according to Regulation (EU) 2017/745 and Council Directive 93/42/EEC respectively:



Manufacturer:

B. Braun Médical | 26, rue Armengaud | 92210 Saint-Cloud | France

The product trademark 'VenaTech' is a registered trademark of B. Braun Melsungen AG or its affiliates.

CE 0459



Distributor:

B. Braun Melsungen AG | Vascular Systems | Sieversufer 8 | 12359 Berlin | Germany

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products.
Reproduction in any form partial or otherwise is not permitted.

6050207 0222/0.5/3